

Centers for Disease Control and Prevention
(CDC)

National Center for Environmental Health
(NCEH)

Division of Laboratory Sciences (DLS)

**NEWBORN SCREENING AND
MOLECULAR BIOLOGY BRANCH
(NSMBB)**

**NEWBORN SCREENING QUALITY
ASSURANCE PROGRAM (NSQAP)
PORTAL**

QUALITY CONTROL USER GUIDE

July 2026

Table of Contents

Portal QC Guide.....2

1. Introduction.....2

2. QC Manual Data Entry.....3

 2.1 QC Manual Data Entry.....3

 2.2 Review Manually Entered Data.....9

 2.3 Submit Manually Entered Data.....16

3. QC Data Upload.....18

 3.1 Pre-filled Template.....22

 3.2 Upload Data.....25

 3.3 Review Uploaded Data.....27

 3.4 Submit Uploaded Data.....30

 3.6 Troubleshooting and FAQs.....31

Portal QC Guide

1. Introduction

The NSMBB Newborn Screening Quality Assurance Program (NSQAP) Portal has been designed to accept participant quality control (QC) data using two options: manual data entry and data upload. Participants are advised to consider their laboratory's internal workflow, processes, and needs, before choosing **one** of the two options for data submission. Additionally, participants are advised to consider the following before deciding on a data entry option:

1. Participants are permitted to only use one data entry option (manual or upload) per event.
2. If the upload option is chosen, data must be consolidated into a single template for upload. Participants are required to download the template, enter data from all QC programs into a single template file and save the file before uploading. If multiple uploads are performed, the previous file will be overwritten.
3. The review feature enables participants to review and edit manually entered or uploaded data in the Portal directly from the review page.

The following pages describe how to enter data into the NSQAP Portal once a data entry option is chosen.


2. QC Manual Data Entry

QC Programs data can be manually entered into the NSQAP Portal from the QC Information section of the Portal.


Note: If you choose to manually enter data, you cannot upload data to the Portal.

2.1 QC Manual Data Entry

1. Select the 'QC Information' button at the top of the page on the toolbar and select the 'QC Information' option.

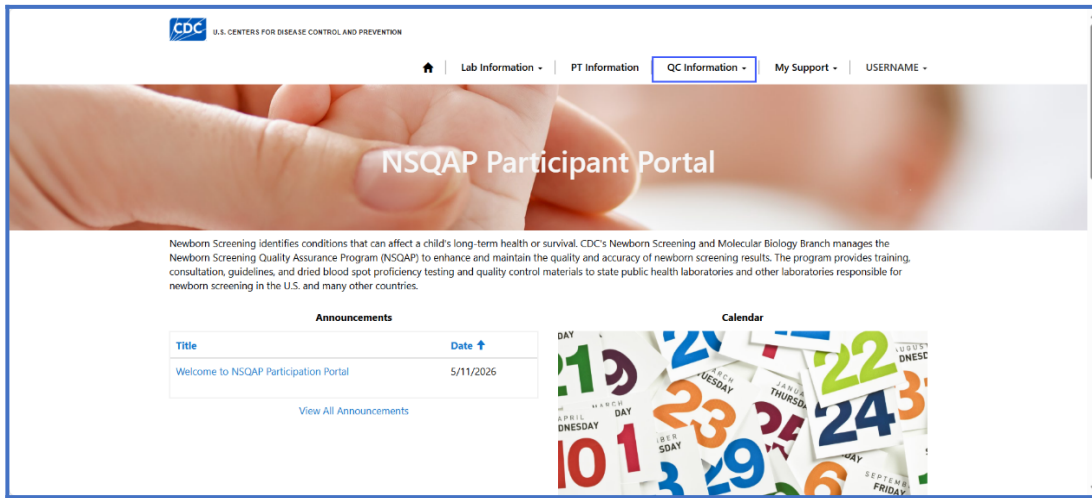

2. Select the 'QC Manual Data Entry' icon on the QC Information homepage.

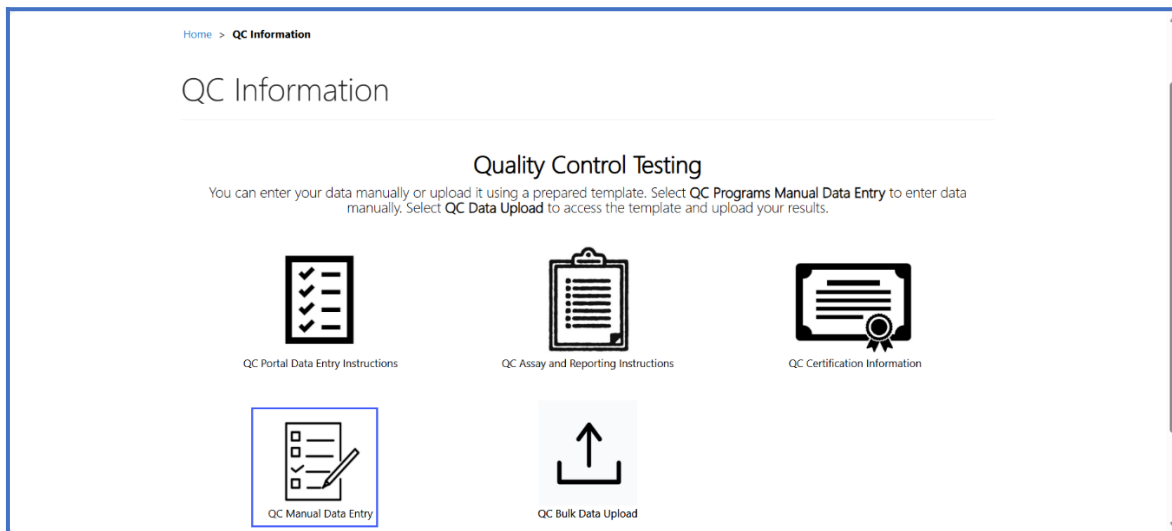

3. Select the 'QC Portal Data Manual Data Entry' icon on the QC Information homepage.



4. Select the QC Program by clicking on the program hyperlink.

QC Programs

Program Name ↑	Created On
17 α-Hydroxyprogesterone + Total Galactose (17OHPQC and TGalQC)	12/10/2023 8:59 PM
2nd Tier Materials (2TM QC)	2/7/2025 5:23 PM
anti-HIV-1 Antibodies in dried blood spots (HIVQC)	12/10/2023 8:59 PM
Duchenne Muscular Dystrophy (DMDQC)	2/1/2025 5:10 AM
Galactose-1-phosphate Uridyltransferase (GALTQC)	12/10/2023 8:59 PM
Immunoreactive Trypsinogen (IRTQC)	12/10/2023 8:59 PM
Lysosomal Storage Disorders (LSDQC)	12/10/2023 8:59 PM
Psychosine (PSYQC)	2/7/2025 5:23 PM

< 1 2 >
*-Required Field.

5. Select the Analyte by clicking on the analyte abbreviation hyperlink.

Analytes List - (DATA ENTRY)

Program Name:
2nd Tier Materials (2TM QC)

Analytes

Abbreviation ↑	Name	Analyte Listing Order	
ALE	Alloisoleucine (ALE)	1	▼
C20-LPC	20:0-Lysophosphatidylcholine (C20-LPC)	2	▼
C22-LPC	22:0-Lysophosphatidylcholine (C22-LPC)	3	▼
C24-LPC	24:0-Lysophosphatidylcholine (C24-LPC)	4	▼
C26-LPC	26:0-Lysophosphatidylcholine (C26-LPC)	5	▼
C3DC	Malonylcarnitine (C3DC)	6	▼
C4OH	Hydroxybutyrylcarnitine (C4OH)	7	▼
CRE	Creatine (CRE)	8	▼
CRN	Creatinine (CRN)	9	▼
EMA	Ethylmalonic Acid (EMA)	10	▼
GUAC	Guanidinoacetic Acid (GUAC)	11	▼
ILE	Isoleucine (ILE)	12	▼
LEU	Leucine (LEU)	13	▼
MA	Malonic Acid (MA)	14	▼
MCA	2-Methylcitric Acid (MCA)	15	▼
MMA	Methylmalonic Acid (MMA)	16	▼
t-HCY	Total Homocysteine (t-HCY)	17	▼
VAL	Valine (VAL)	18	▼

6. Select the 'Magnifying Glass' icon on the 'Method' field.

QC Data Entry

[View Summary](#)

[SAVE METHOD](#)

General

Analyte Name
Alloisoleucine (ALE)

Method *



Runs

7. A new window will appear listing methods for the analyte. Click on the method, and the row will highlight with a check mark on the left side. Click the **'Select'** button at the bottom of the window to choose the method.

Lookup records

Search

Choose one record and click Select to continue

✓

Method Name ↑

☐	LC-MS/MS negative ion mode
☐	LC-MS/MS non-kit
☐	LC-MS/MS positive ion mode
☐	Other
☐	Ultra Performance Liquid Chromatography

Select

Cancel

Remove value

8. Select the **'Save Method'** button to save the method for the analyte.

QC Data Entry

View Summary

SAVE METHOD

General


Analyte Name

Alloisoleucine (ALE)

Method *

Runs

- Select the 'OK' button when prompted "are you sure you want to save the method?".

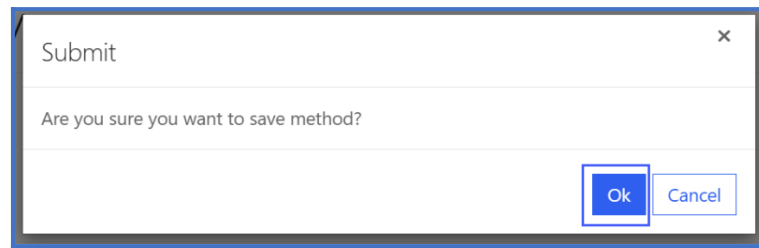

- To add QC data, select the run hyperlink in the 'Runs' summary table.

Runs

The summary table will display data with two decimal places, regardless of the values entered. Results entered as <LOD will not be displayed in the summary table.

[Download Summary](#)

	Analyte	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot_C	Replicate_1C	Replicate_2C	Lot_D	Re
1	Alloisoleucine (ALE)	1111111			222222			333333			444444	
2	Alloisoleucine (ALE)	1111111			222222			333333			444444	
3	Alloisoleucine (ALE)	1111111			222222			333333			444444	
4	Alloisoleucine (ALE)	1111111			222222			333333			444444	
5	Alloisoleucine (ALE)	1111111			222222			333333			444444	

11. A data entry window will appear. Enter a quantitative value or select the <LOD checkbox for both replicates per lot then select the **'SAVE RUN DATA'** button.

View details

Analyte: Alloisoleucine (ALE)

1

SAVE RUN DATA

Report data in µmol/L

Data cannot not be saved until all replicates have been completed

Values

1111111

Replicate 1A *

☐ <LOD

Replicate 2A *

☐ <LOD

NOTE: To save run data, both replicates must have a quantitative value or <LOD selected.

12. Data will appear in the summary table when the run is saved.

Runs

The summary table will display data with two decimal places, regardless of the values entered. Results entered as <LOD will not be displayed in the summary table.

Download Summary

RUNS	Analyte	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot_C	Replicate_1C	Replicate_2C	Lot_D	Re
1	Alloisoleucine (ALE)	1111111	1.00	1.00	222222	1.00	1.00	333333	1.00	1.00	444444	1.00
2	Alloisoleucine (ALE)	1111111			222222			333333			444444	
3	Alloisoleucine (ALE)	1111111			222222			333333			444444	
4	Alloisoleucine (ALE)	1111111			222222			333333			444444	
5	Alloisoleucine (ALE)	1111111			222222			333333			444444	

13. Analyte specific data can be entered, reviewed, and edited on this page.

NOTE: Export results to an Excel spreadsheet by selecting the 'Download Summary' button.
(Optional)

Runs

The summary table will display data with two decimal places, regardless of the values entered. Results entered as <LOD will not be displayed in the summary table.


[Download Summary](#)

RUNS ↑	Analyte	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot_C	Replicate_1C	Replicate_2C	Lot_D	Re
1	Alloisoleucine (ALE)	1111111	1.00	1.00	222222	1.00	1.00	333333	1.00	1.00	444444	1.0
2	Alloisoleucine (ALE)	1111111			222222			333333			444444	
3	Alloisoleucine (ALE)	1111111			222222			333333			444444	
4	Alloisoleucine (ALE)	1111111			222222			333333			444444	
5	Alloisoleucine (ALE)	1111111			222222			333333			444444	

2.2 Review Manually Entered Data

All manually entered QC Programs data can be reviewed by accessing the 'QC Manual Review' page on the QC Information page.

1. Select the 'QC Manual Data Entry' icon from the QC Information page.

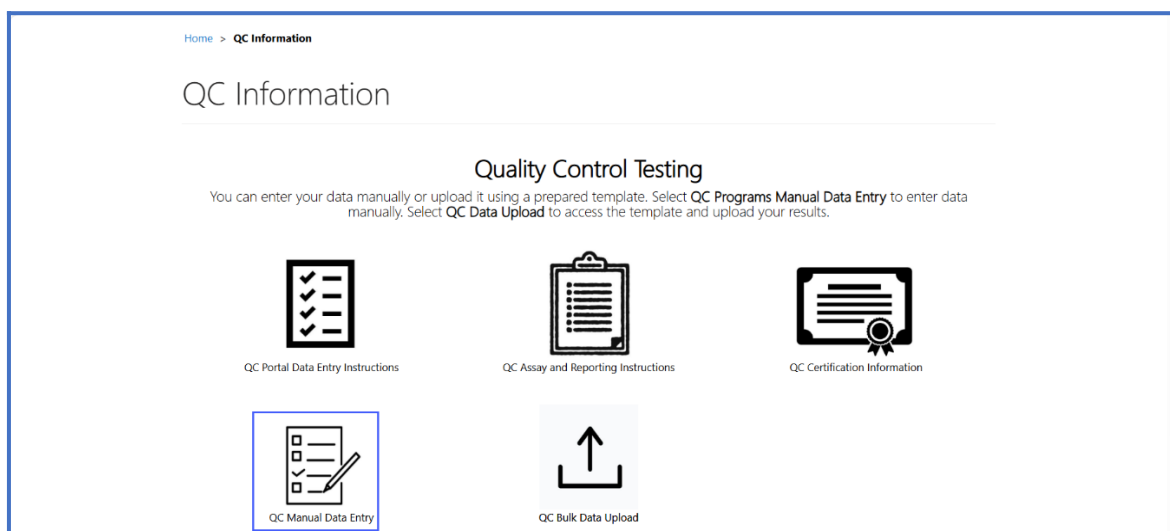

2. Select the 'QC Manual Review' icon from the QC Information page.

Quality Control Testing-Manual Entry

You can enter your data manually. Select **QC Programs Manual Data Entry** to enter data manually



QC Portal Data Manual
Data Entry



QC Manual Review



QC Manual Submit

3. On the QC Manual Review page, results can be reviewed and edited.

QC Manual Review

To filter the data by program, select the check box(es) next to the program and click the "Apply" button in the table below.

Filter by Program Name

☐ 17 α -Hydroxyprogesterone + Total Galactose (17OHPQC and TGalQC)

☒ 2nd Tier Materials (2TMQC)

☐ Duchenne Muscular Dystrophy (DMDQC)

☐ Galactose-1-phosphate Uridyltransferase (GALTQC)

☐ Immunoreactive Trypsinogen (IRTQC)

☐ Lysosomal Storage Disorders (LSDQC)

☐ Psychosine (PSYQC)

☐ Second-tier Congenital Adrenal Hyperplasia by LC-MS/MS (CAHQC)

☐ Second-tier Maple Syrup Urine Disease and Phenylketonuria by LC-MS/MS (MSUD-PKUQC)

☐ Second-tier Methylmalonic /Propionic Acidemia and Homocystinuria by LC-MS/MS (MMA-HCYQC)

☐ Tandem MS 1 (MSMS1QC)

☐ Thyroid-Stimulating Hormone (TSHQC)

☐ Thyroxine (T4QC)

[Less](#) ▼

[Apply](#)

[Download](#)

Run ↑	Method	Method Code	Analyte	Analyte_code	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot
1			PSY	501	PSY	111			222			333
1			Creatine (CRE)	509	CRE	1111111			222222			333
1			CK-MM	500	CK-MM	111			222			333
1			Creatinine (CRN)	510	CRN	1111111			222222			333

- To review QC data, navigate through the results table using the navigation features located on the right side and the bottom of the table.

[Download](#)

Run ↑	Method	Method Code	Analyte	Analyte_code	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot
1			PSY	501	PSY	111			222			333
1			Creatine (CRE)	509	CRE	1111111			222222			333
1			CK-MM	500	CK-MM	111			222			333
1			Creatinine (CRN)	510	CRN	1111111			222222			333
1			Malonylcarnitine (C3DC)	507	C3DC	1111111			222222			333
1	LC-MS/MS negative ion mode		240- lysophosphatidylcholine (C24-LPC)	505	C24- LPC	1111111			222222			333
1	Derivatized - MS/MS ClenSpot® Complete Kit RECIFE		200- lysophosphatidylcholine (C20-LPC)	503	C20- LPC	1111111			222222			333
1			220- lysophosphatidylcholine (C22-LPC)	504	C22- LPC	1111111			222222			333
1			240- lysophosphatidylcholine (C24-LPC)	505	C24- LPC	1111111			222222			333
1	Derivatized - MS/MS NeoGram PerkinElmer		200- lysophosphatidylcholine (C20-LPC)	503	C20- LPC	1111111			222222			333
1			260- lysophosphatidylcholine (C26-LPC)	506	C26- LPC	1111111			222222			333
1			CK-MM	500	CK-MM	111			222			333
1	LC-MS/MS non-ikt		Allisoleucine (ALE)	502	ALE	1111111	1.00	1.00	222222	1.00	1.00	333
1			220- lysophosphatidylcholine (C22-LPC)	504	C22- LPC	1111111			222222			333
1	Other		PSY	501	PSY	111			222			333

<
1
2
3
4
5
6
7
8
...
66
>

Form Approved
OMB No. 0502-1389
Exp. Date 01/31/2026

- To edit data, select the hyperlinked run number of the row of interest.

Run ↑	Method	Method Code	Analyte	Analyte_code	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot
1			PSY	501	PSY	111			222			333
1			Creatine (CRE)	509	CRE	1111111			222222			333
1			CK-MM	500	CK-MM	111			222			333
1			Creatinine (CRN)	510	CRN	1111111			222222			333

6. The analyte data entry page will appear. The data can be edited and saved.

View details
 ×

Analyte: Alloisoleucine (ALE)

1

SAVE RUN DATA

Report data in $\mu\text{mol/L}$

Data cannot not be saved until all replicates have been completed

Values

1111111

Replicate 1A *

☐ <LOD

Replicate 2A *

☐ <LOD

7. After selecting 'Save Run Data' you will be directed back to the QC Review page to continue reviewing your data.

QC Manual Review

To filter the data by program, select the check box(es) next to the program and click the "Apply" button in the table below.

Filter by Program Name

☐ 17 α -Hydroxyprogesterone + Total Galactose (17OHPQC and TGalQC)

☐ 2nd Tier Materials (2TMQC)

☐ Duchenne Muscular Dystrophy (DMDQC)

☐ Galactose-1-phosphate Uridyltransferase (GALTQC)

☐ Immunoreactive Trypsinogen (IRTQC)

☐ Lysosomal Storage Disorders (LSDQC)

☐ Psychosine (PSYQC)

☐ Second-tier Congenital Adrenal Hyperplasia by LC-MS/MS (CAHQC)

☐ Second-tier Maple Syrup Urine Disease and Phenylketonuria by LC-MS/MS (MSUD-PKUQC)

☐ Second-tier Methylmalonic /Propionic Acidemia and Homocystinuria by LC-MS/MS (MMA-HCYQC)

☐ Tandem MS 1 (MSMS1QC)

☐ Thyroid-Stimulating Hormone (TSHQC)

☐ Thyroxine (T4QC)

[Less](#) ▼

[Apply](#)

[Download](#)

Run ↑	Method	Method Code	Analyte	Analyte_code	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot
1			PSY	501	PSY	111			222			333
1			Creatine (CRE)	509	CRE	1111111			222222			333
1			CK-MM	500	CK-MM	111			222			333
1			Creatinine (CRN)	510	CRN	1111111			222222			333

8. Review QC data offline by selecting 'Download' to generate an Excel spreadsheet of the data entered.

QC Manual Review

To filter the data by program, select the check box(es) next to the program and click the "Apply" button in the table below.

Filter by Program Name

☐ 17 α -Hydroxyprogesterone + Total Galactose (17OHPQC and TGalQC)

☐ 2nd Tier Materials (2TMQC)

☐ Duchenne Muscular Dystrophy (DMDQC)

☐ Galactose-1-phosphate Uridyltransferase (GALTQC)

☐ Immunoreactive Trypsinogen (IRTQC)

☐ Lysosomal Storage Disorders (LSDQC)

☐ Psychosine (PSYQC)

☐ Second-tier Congenital Adrenal Hyperplasia by LC-MS/MS (CAHQC)

☐ Second-tier Maple Syrup Urine Disease and Phenylketonuria by LC-MS/MS (MSUD-PKUQC)

☐ Second-tier Methylmalonic /Propionic Acidemia and Homocystinuria by LC-MS/MS (MMA-HCYQC)

☐ Tandem MS 1 (MSMS1QC)

☐ Thyroid-Stimulating Hormone (TSHQC)

☐ Thyroxine (T4QC)

[Less](#)

[Apply](#)

[Download](#)

Run	Method	Method Code	Analyte	Analyte_code	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot
1			PSY	501	PSY	111			222			333
1			Creatine (CRE)	509	CRE	1111111			222222			333
1			CK-MM	500	CK-MM	111			222			333
1			Creatinine (CRN)	510	CRN	1111111			222222			333

9. QC data can be filtered by program. Check the box(es) next to all relevant program names then select 'Apply'.

QC Manual Review

To filter the data by program, select the check box(es) next to the program and click the "Apply" button in the table below.

Filter by Program Name

☐ 17 α -Hydroxyprogesterone + Total Galactose (17OHPQC and TGalQC)
 ☒ 2nd Tier Materials (2TMQC)
 ☐ Duchenne Muscular Dystrophy (DMDQC)
 ☐ Galactose-1-phosphate Uridyltransferase (GALTQC)
 ☐ Immunoreactive Trypsinogen (IRTQC)
 ☐ Lysosomal Storage Disorders (LSDQC)
 ☐ Psychosine (PSYQC)
 ☐ Second-tier Congenital Adrenal Hyperplasia by LC-MS/MS (CAHQC)
 ☐ Second-tier Maple Syrup Urine Disease and Phenylketonuria by LC-MS/MS (MSUD-PKUQC)
 ☐ Second-tier Methylmalonic /Propionic Acidemia and Homocystinuria by LC-MS/MS (MMA-HCYQC)
 ☐ Tandem MS 1 (MSMS1QC)
 ☐ Thyroid-Stimulating Hormone (TSHQC)
 ☐ Thyroxine (T4QC)

[Less](#)


Apply

Download


Run	Method	Method Code	Analyte	Analyte_code	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot
1		PSY		501	PSY	111			222			333
1			Creatine (CRE)	509	CRE	1111111			222222			333
1		CK-MM		500	CK-MM	111			222			333
1			Creatinine (CRN)	510	CRN	1111111			222222			333

2.3 Submit Manually Entered Data


1. After reviewing data, navigate to the '**QC Manual Data Entry**' page on the '**QC Information**' page to submit all data at once.

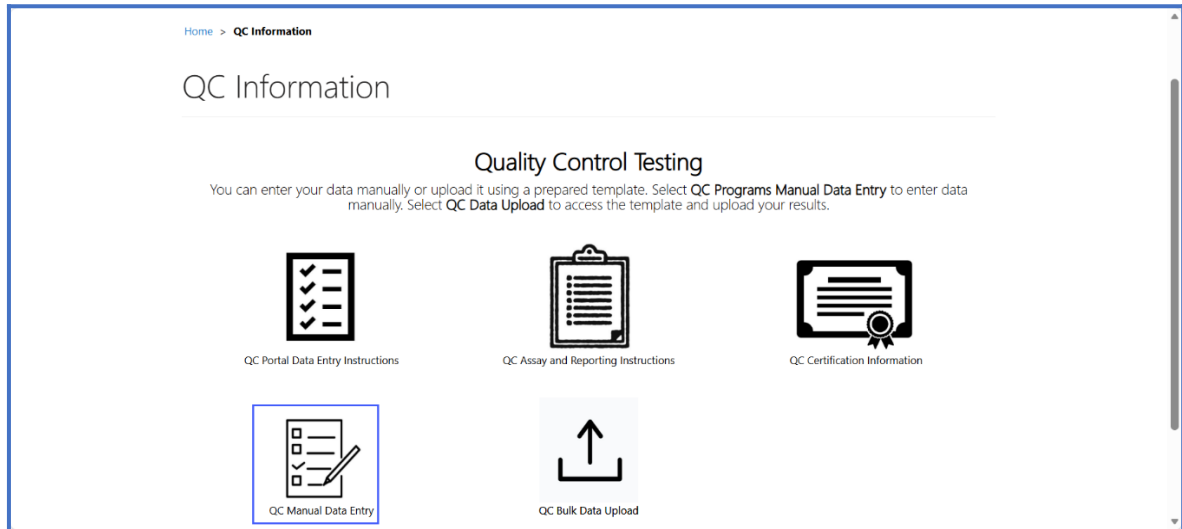

2. Click the '**QC Manual Submit**' button on the '**Quality Control Testing-Manual Entry**' page to submit all data at once.

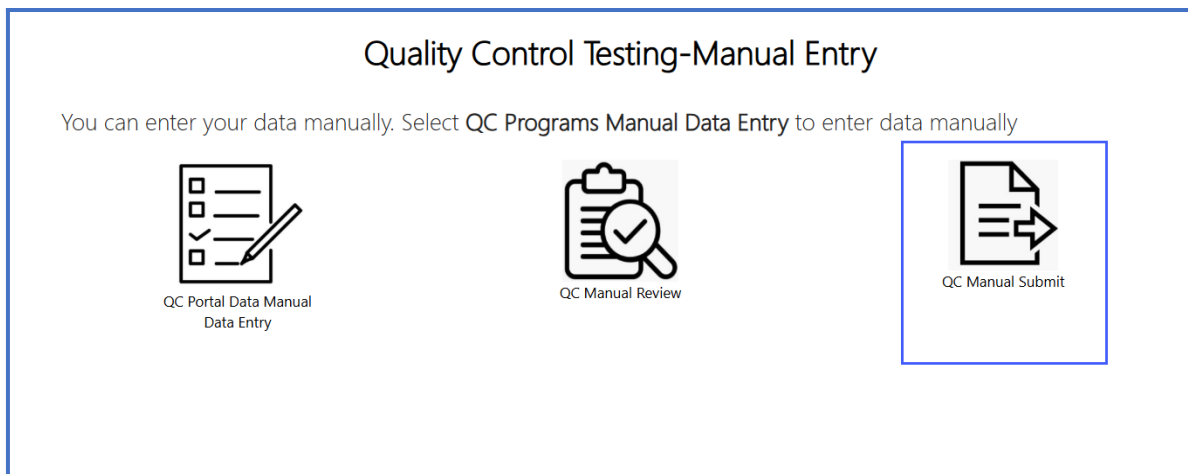

3. Select 'Submit'.

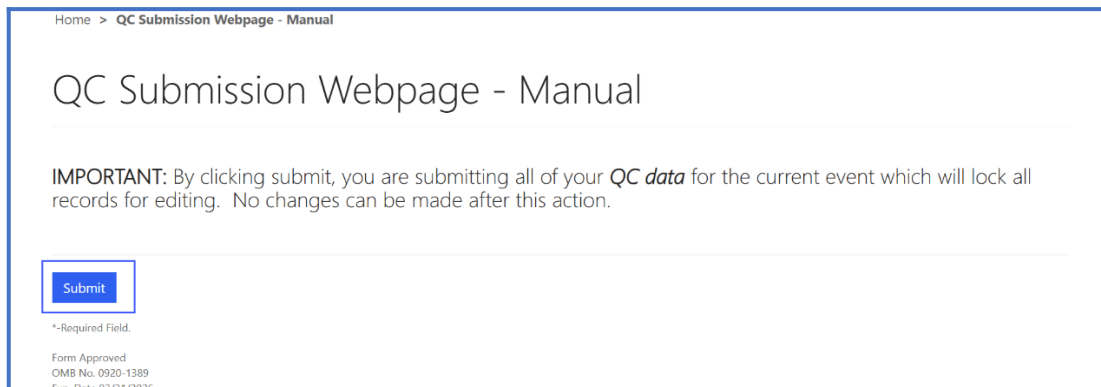

Home > QC Submission Webpage - Manual

QC Submission Webpage - Manual


IMPORTANT: By clicking submit, you are submitting all of your **QC data** for the current event which will lock all records for editing. No changes can be made after this action.

[Submit](#)

*-Required Field.

Form Approved
OMB No. 0920-1389

4. A confirmation page will appear once the submission is complete.

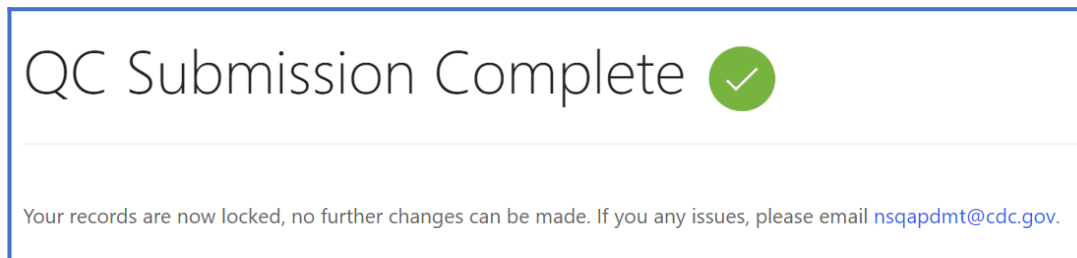

QC Submission Complete


Your records are now locked, no further changes can be made. If you any issues, please email nsqapdmt@cdc.gov.

3. QC Data Upload


Data for QC Programs can be uploaded to the NSQAP Portal by accessing the 'QC Data Upload' area.

Note: If you choose to upload data, you cannot manually enter data into the Portal.

1. Select the '**QC Bulk Data Upload**' icon on the QC information page.

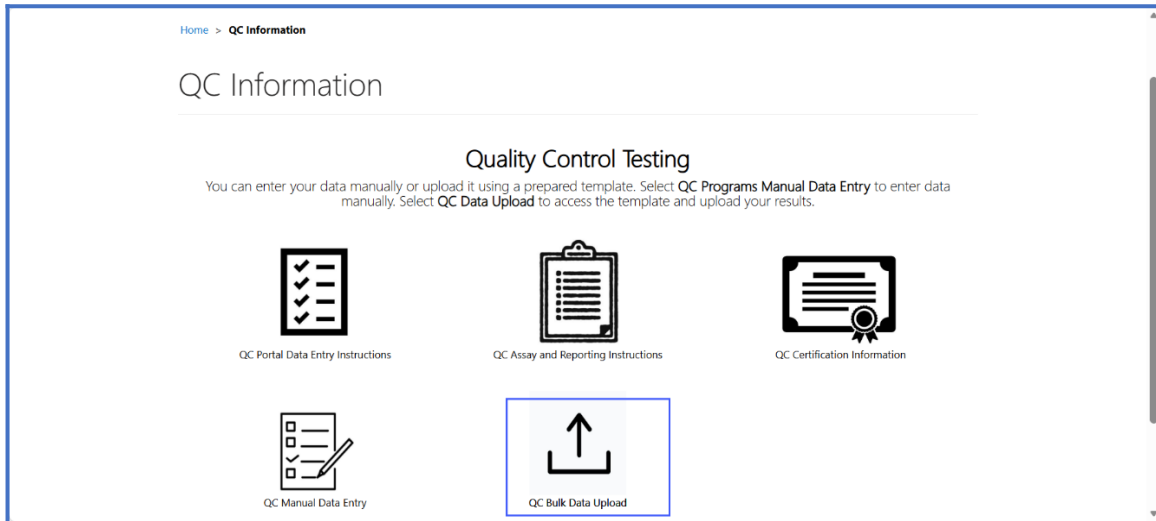

2. The QC Data Upload Information homepage contains several icons to navigate to the various QC upload sections.

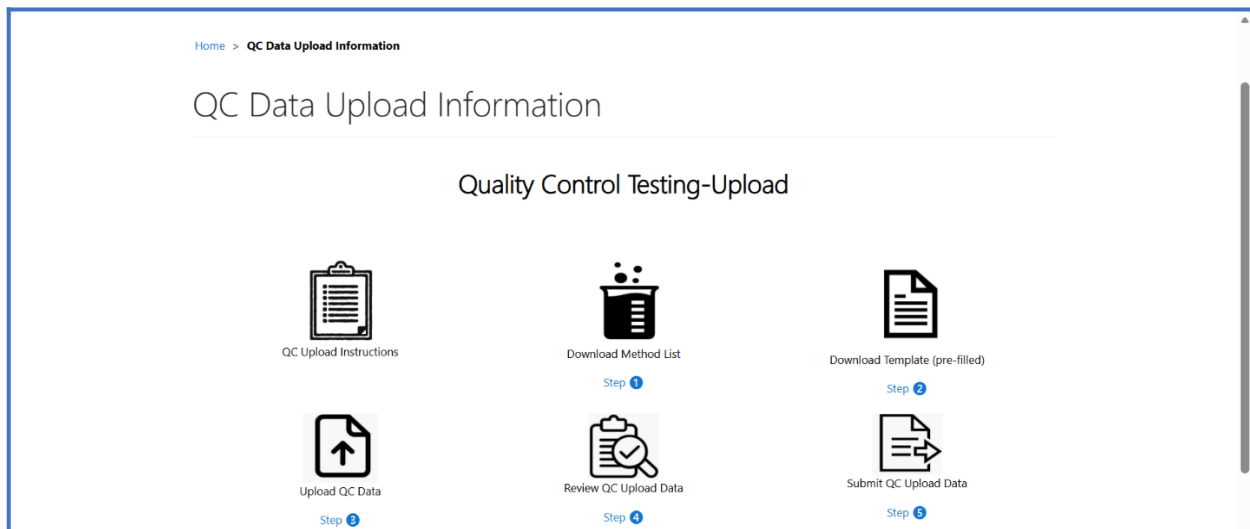

3. QC Upload Instructions – Instructions for data upload process.

Upload Instructions

Steps for Uploading QC Template:

1. Download the Pre-filled template option.
2. Download and **READ** the QC Upload Instructions.
3. Download the Method List.
4. Enter all QC Programs data into a single consolidated template file. Report exactly 5 runs (10 data points) for each analyte. This is required for the template to be accepted.
 - When using the pre-filled template, provide the following information: method code, other method name (if applicable), and results.
5. Do not alter column names, column order, or format of your template file. Any changes to the structure or format of the template will inactivate the template and the upload will fail.
6. Leave blank any inapplicable columns or fields. Add <LOD for results when necessary.
7. Upload the completed template for **ALL PROGRAMS** you are reporting. No partial results will be accepted at this time
8. **If you upload more than once, your previous upload will be overwritten.**
9. Use the 'Review QC Upload Data' page to review and edit entries after uploading the completed template.
10. If you have trouble uploading your completed template please email NSQAPDMT@cdc.gov for assistance.

4. Download Method List – List of methods and their associated codes. (Method codes are needed for data upload.)

Method List

- **DOWNLOAD AND READ INSTRUCTIONS ON HOW TO USE THE QC UPLOAD TEMPLATE.**

- Use the method table below as a reference for entering method codes when using the pre-filled and blank QC templates.
- Method codes can be found by navigating through the pages of the table or searching by searching using the search bar.
- Note: Use an asterisk to search method names and codes using partial text (e.g. Enter *neobase to search for Non-derivatized - MS/MS NeoBase™ PerkinElmer).



 Download

Method Name ↑

Method Code

50hr Reagent Kit Spotcheck® BIOT Astoria-Pacific, ERU (1µmol/dL/90min)

360

50hr Reagent Kit Spotcheck® G6PD Astoria-Pacific, µmol/L blood

320

50hr Reagent Kit Spotcheck® GALT Astoria-Pacific, µmol/L blood

240

50hr Reagent Kit Spotcheck® TGal Astoria-Pacific

200

Accuwell™ EIA TSH Neo-Genesis

120

AutoDELFIA® Neonatal 17OHP PerkinElmer

160

AutoDELFIA® Neonatal hTSH PerkinElmer

121

AutoDELFIA® Neonatal IRT PerkinElmer

280

5. Download Template (pre-filled) – Downloadable template that is pre-filled with QC Lot Numbers, analytes, and analyte codes.

Download QC Template - (Pre-filled)

Note: You will see only the program(s) in which you are enrolled.

To obtain a pre-filled template click download (obtain a current template for each QC event). A completed template can be replaced up until the data is Submitted. The most recently uploaded file will be displayed. Data cannot be edited or replaced after the file is Submitted.

Do not change any of the column names as your data will not be uploaded correctly and the file will not be accepted. For more information please go to the resource area and review/download directions.

 Download

RUN	Method	OTHER_METHOD	Analyte_code	Analyte ↓	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot_C
	Code											
1			64	Total Galactose (TGal)	TGal	A0101			A0102			A0103
2			64	Total Galactose (TGal)	TGal	A0101			A0102			A0103
3			64	Total Galactose (TGal)	TGal	A0101			A0102			A0103

6. Upload QC Data – Page where completed templates are uploaded.

QC Data Upload

Steps for Uploading QC Template:

1. Download one of the template options "Blank" or "Pre-filled".
2. Download and **READ** the QC Uploading Procedure.
3. Download the Method List.
4. Download the Analyte List if using "Blank" Template.
5. Enter all QC Programs data into a single consolidated template file. Report exactly 5 runs (10 data points) for each analyte. This is required for the template to be accepted.
 - For the blank template, provide the following information: method code, other method name (if applicable), analyte code, run number (1,2,3,4,5), lot number, and results.
 - For the pre-filled template, provide the following information: method code, other method name (if applicable), and results.
6. Do not alter column names, column order, or format of your template file. Any changes to the structure or format of the template will inactivate the template and the upload will fail.
7. Leave blank any inapplicable columns or fields. Add <LOD for results when necessary.
8. Upload the completed template for **ALL QC PROGRAMS** you are reporting. No partial results will be accepted at this time.
9. **If you upload more than once, your previous file will be overwritten.**
10. Use the 'Review QC Upload Data' page to review and edit entries (if necessary) after uploading the completed template.
11. If you have trouble uploading your completed template please email NSQAPDMT@cdc.gov for assistance.

Attach the Excel Template file *

No file chosen

7. Review QC Upload Data – Review uploaded data before submission.

QC Review and Edit Upload Data

To filter the data by program, select the check box(es) next to the program and click the "Apply" button in the table below.

Filter by Program Name

- ☐ 17 α -Hydroxyprogesterone + Total Galactose (17OHPQC and TGalQC)
- ☐ Galactose-1-phosphate Uridyltransferase (GALTQC)
- ☐ Immunoreactive Trypsinogen (IRTQC)
- ☐ Lysosomal Storage Disorders (LSDQC)
- ☐ Second-tier Congenital Adrenal Hyperplasia by LC-MS/MS (CAHQC)
- ☐ Second-tier Maple Syrup Urine Disease and Phenylketonuria by LC-MS/MS (MSUD-PKUQC)

More ▼

Apply

Download

RUN ↑	Method	Method Code	OTHER_METHOD	Analyte	Analyte Code	Analyte Abbreviation	Lot_A	Replicate_1A	Replicate_2A
1	Colormetric BIOT non- kit, qualitative	365		α -L- Iduronidase (IDUA)	88	IDUA	A1808	366.00	359.00
2	Colormetric BIOT non- kit, qualitative	365		α -L- Iduronidase (IDUA)	88	IDUA	A1808	366.00	367.00

8. Submit QC Upload Data – Page to submit uploaded data.

Home > QC Submission Webpage - Upload


QC Submission Webpage - Upload

IMPORTANT: By clicking submit, you are submitting all of your *QC data* for the current event which will lock all records for editing. No changes can be made after this action.

Submit

3.1 Pre-filled Template

1. Select the **'Download Template (pre-filled)'** icon from the QC Data Upload Information page.

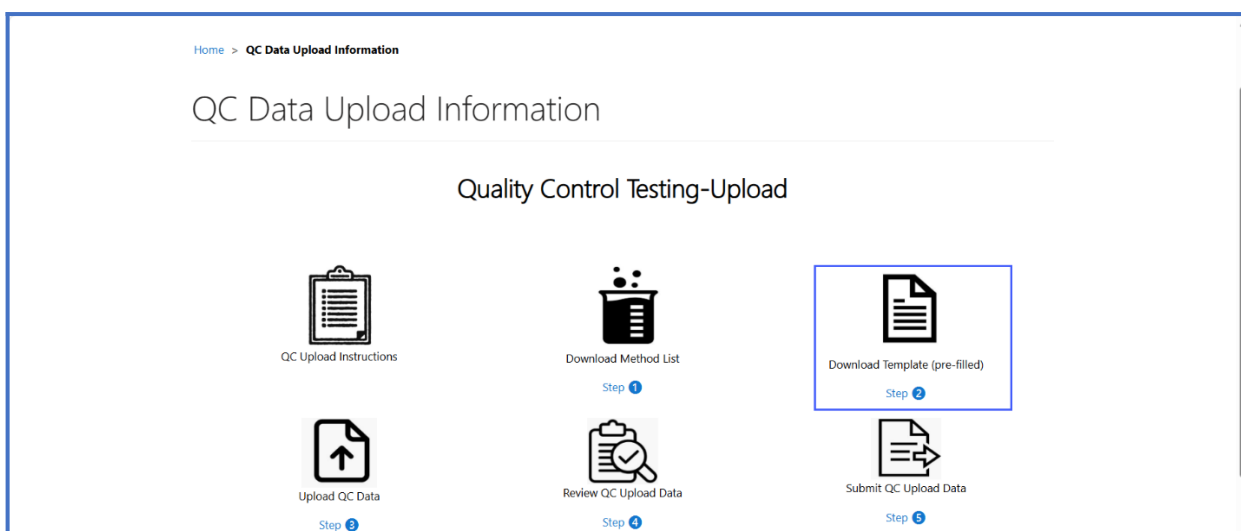

2. Select the **'Download'** button from the Pre-filled Template Page.

Download QC Template - (Pre-filled)


Note: You will see only the program(s) in which you are enrolled.

To obtain a pre-filled template click download (obtain a current template for each QC event). A completed template can be replaced up until the data is Submitted. The most recently uploaded file will be displayed. Data cannot be edited or replaced after the file is Submitted. Do not change any of the column names as your data will not be uploaded correctly and the file will not be accepted. For more information please go to the resource area and review/download directions.

[Download](#)

RUN	Method Code	OTHER_METHOD	Analyte_code	Analyte ↓	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot_C
1			64	Total Galactose (TGal)	TGal	A0101			A0102			A0103
2			64	Total Galactose (TGal)	TGal	A0101			A0102			A0103
3			64	Total Galactose (TGal)	TGal	A0101			A0102			A0103

3. Select the 'Open' button on the downloaded Excel file at the bottom of the page.

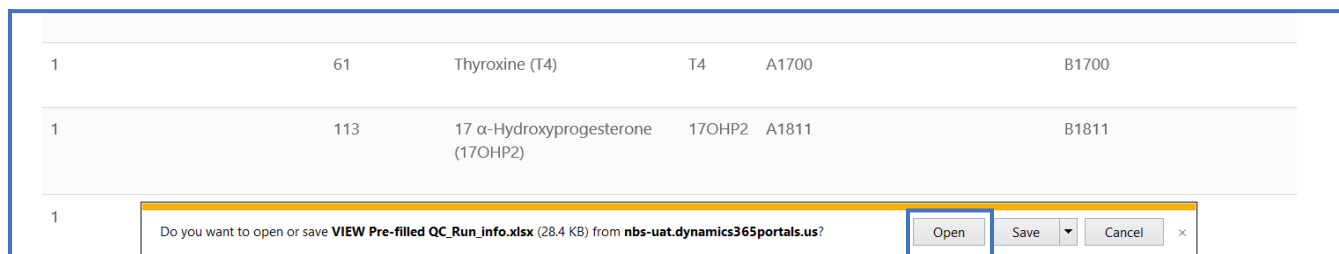

Note: The Excel download may look different depending on the browser used. This screenshot is from the Internet Explorer browser.

4. Complete the Excel template with data for submission and save the file.

To complete the Pre-filled template, you will need to include:

1. Method Code
2. Replicate Values (Insert "<LOD" where necessary)

	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P
1	RUN	Method Code	OTHER_METHOD	Analyte_code	Analyte	Abbr	Lot_A	Replicate_1A	Replicate_2A	Lot_B	Replicate_1B	Replicate_2B	Lot_C	Replicate_1C	Replicate_2C	Lot_D
2	1	403	103	103	Tiglylcarnitine (C5:1	A1815	0.13	0.33 B1815		0.45	0.43 C1815		1.35	1.32 D1815		
3	2	403	103	103	Tiglylcarnitine (C5:1	A1815	0.12	0.12 B1815		0.49	0.43 C1815		1.24	1.34 D1815		
4	3	403	103	103	Tiglylcarnitine (C5:1	A1815	0.11	0.09 B1815		0.45	0.46 C1815		1.33	1.43 D1815		
5	4	403	103	103	Tiglylcarnitine (C5:1	A1815	0.11	0.1 B1815		0.42	0.47 C1815		1.39	1.32 D1815		
6	5	403	103	103	Tiglylcarnitine (C5:1	A1815	0.12	0.11 B1815		0.47	0.51 C1815		1.43	1.42 D1815		
7	1	403	104	104	Tetradecenoylc C14:1	A1815	0.1	0.07 B1815		0.49	0.48 C1815		1.45	1.38 D1815		
8	2	403	104	104	Tetradecenoylc C14:1	A1815	0.06	0.07 B1815		0.49	0.53 C1815		1.38	1.48 D1815		
9	3	403	104	104	Tetradecenoylc C14:1	A1815	0.08	0.08 B1815		0.51	0.51 C1815		1.45	1.47 D1815		
10	4	403	104	104	Tetradecenoylc C14:1	A1815	0.09	0.07 B1815		0.51	0.56 C1815		1.53	1.42 D1815		
11	5	403	104	104	Tetradecenoylc C14:1	A1815	0.09	0.08 B1815		0.45	0.49 C1815		1.49	1.35 D1815		
12	1	403	105	105	Creatine (CRE) CRE	A1815	185.57	137.51 B1815		275.83	272.05 C1815		469.75	462.38 D1815		
13	2	403	105	105	Creatine (CRE) CRE	A1815	204.8	208.6 B1815		330.9	301.5 C1815		479.5	487.9 D1815		
14	3	403	105	105	Creatine (CRE) CRE	A1815	186.05	190.15 B1815		290.68	301.17 C1815		473.97	467.75 D1815		
15	4	403	105	105	Creatine (CRE) CRE	A1815	200.76	188.69 B1815		299.16	297.61 C1815		473.3	483.02 D1815		
16	5	403	105	105	Creatine (CRE) CRE	A1815	222.14	208.44 B1815		296.15	322.09 C1815		505.45	474.55 D1815		
17	1	403	106	106	Guanidinoacetyl GUAC	A1815	1.34	2.64 B1815		5.87	5.97 C1815		11.3	10.93 D1815		
18	2	403	106	106	Guanidinoacetyl GUAC	A1815	1.62	1.58 B1815		7.21	6.61 C1815		11.87	11.95 D1815		
19	3	403	106	106	Guanidinoacetyl GUAC	A1815	1.19	1.32 B1815		6.16	6.3 C1815		10.91	11.18 D1815		
20	4	403	106	106	Guanidinoacetyl GUAC	A1815	1.31	1.28 B1815		6.13	6.45 C1815		11.7	12.01 D1815		
21	5	403	106	106	Guanidinoacetyl GUAC	A1815	1.58	1.49 B1815		6.45	6.66 C1815		12.18	11.38 D1815		
22	1	403	107	107	Creatinine (CRN) CRN	A1815	43.06	40.14 B1815		74.92	74.5 C1815		124.48	120.13 D1815		
23	2	403	107	107	Creatinine (CRN) CRN	A1815	39.75	41.19 B1815		81.65	71.52 C1815		116.7	121.98 D1815		
24	3	403	107	107	Creatinine (CRN) CRN	A1815	39.61	40.28 B1815		77.74	82.25 C1815		122.72	124.27 D1815		
25	4	403	107	107	Creatinine (CRN) CRN	A1815	42.21	40.76 B1815		81.65	83.44 C1815		127.02	127.05 D1815		
26	5	403	107	107	Creatinine (CRN) CRN	A1815	43.54	40.34 B1815		71.83	79.7 C1815		128.78	126.39 D1815		
27	1	408	108	108	20:0-Lysophosp C20-LPC	A1815	0.38	0.35 B1815		0.59	0.55 C1815		0.94	0.9 D1815		
28	2	408	108	108	20:0-Lysophosp C20-LPC	A1815	0.34	0.35 B1815		0.54	0.53 C1815		0.94	0.95 D1815		
29	3	408	108	108	20:0-Lysophosp C20-LPC	A1815	0.34	0.31 B1815		0.56	0.55 C1815		1.07	1.05 D1815		

The table below explains each of the columns in the template.

Column name	Description	Values	Example Data
RUN	Shows the run number	Whole number (INTEGER) from 1-5	1
Method Code	Code for the method used	See Method list webpage	360
OTHER_METHOD	Name of method not in list	Type name of method	MYNewMethod
Analyte_code	Code for the analyte	See Analyte list webpage	203
Lot_A	Lot number for Lot A	Name of lot provided for each event	A1811
Replicate_1A	Result for replicate 1 of Lot A	Decimal numbers with periods or <LOD	1.12
Replicate_2A	Result for replicate 2 of Lot A	Decimal numbers with periods or <LOD	1.12
Lot_B	Lot number for Lot B	Name of lot provided each event	B1811
Replicate_1B	Result for replicate 1 of Lot B	Decimal numbers with periods or <LOD	1.12
Replicate_2B	Result for replicate 2 of Lot B	Decimal numbers with periods or <LOD	1.12
Lot_C	Lot number for Lot C	Name of lot provided each event	C1811
Replicate_1C	Result for replicate 1 of Lot C	Decimal numbers with periods or <LOD	1.12
Replicate_2C	Replicate 2 of Lot C	Decimal numbers with periods or <LOD	1.12
Lot_D	Lot number for Lot D	Name of lot provided each event	D1811
Replicate_1D	Result for replicate 1 of Lot D	Decimal numbers with periods or <LOD	1.12
Replicate_2D	Result for replicate 2 of Lot D	Decimal numbers with periods or <LOD	1.12
Lot_E	Lot number for Lot E	Name of lot provided each event	E1811
Replicate_1E	Result for replicate 1 of Lot E	Decimal numbers with periods or <LOD	1.12
Replicate_2E	Result for replicate 2 of Lot E	Decimal numbers with periods or <LOD	1.12


Note: If you choose to complete the Pre-filled Template, you cannot complete the Blank Template and submit it.

Note: Method Codes and Analyte Codes are found on the 'Download Method List' and the 'Download Analyte List' pages, respectively.

Note: Do not change any of the column headers. Headers must be exact for information to be uploaded properly.

3.2 Upload Data

1. Navigate back to the 'QC Data Upload Information Page' and select the 'Upload QC Data' icon.

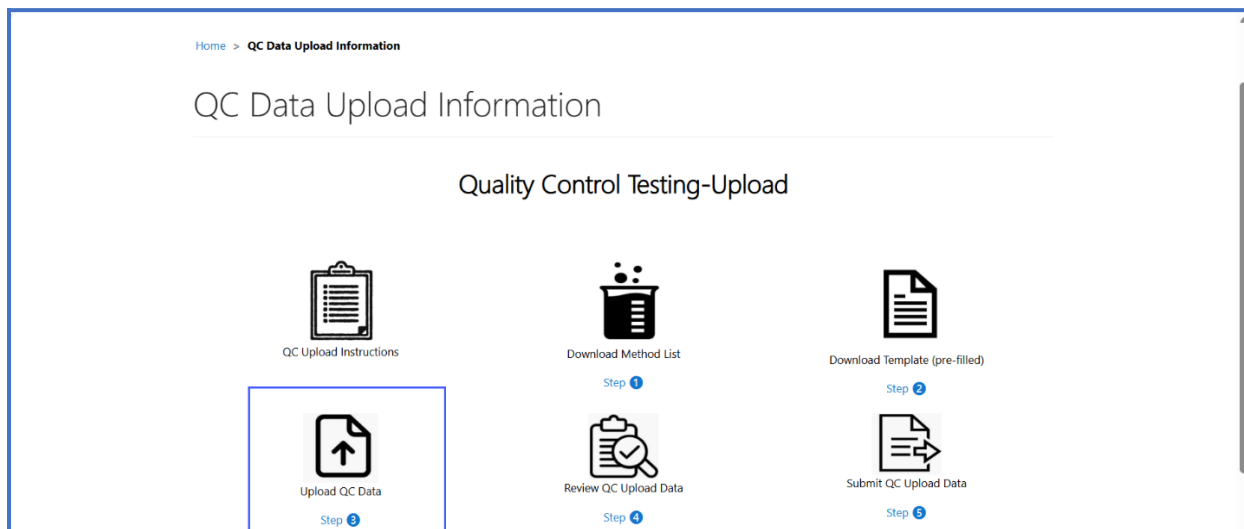

2. Select the 'Choose File' button on the 'QC Data Upload' page and select the completed template for upload.

QC Data Upload

Steps for Uploading QC Template:

1. Download one of the template options "Blank" or "Pre-filled".
2. Download and **READ** the QC Uploading Procedure.
3. Download the Method List.
4. Download the Analyte List if using "Blank" Template.
5. Enter all QC Programs data into a single consolidated template file. Report exactly 5 runs (10 data points) for each analyte. This is required for the template to be accepted.
 - For the blank template, provide the following information: method code, other method name (if applicable), analyte code, run number (1,2,3,4,5), lot number, and results.
 - For the pre-filled template, provide the following information: method code, other method name (if applicable), and results.
6. Do not alter column names, column order, or format of your template file. Any changes to the structure or format of the template will inactivate the template and the upload will fail.
7. Leave blank any inapplicable columns or fields. Add <LOD for results when necessary.
8. Upload the completed template for **ALL QC PROGRAMS** you are reporting. No partial results will be accepted at this time.
9. **If you upload more than once, your previous file will be overwritten.**
10. Use the 'Review QC Upload Data' page to review and edit entries (if necessary) after uploading the completed template.
11. If you have trouble uploading your completed template please email NSQAPDMT@cdc.gov for assistance.

Attach the Excel Template file *

No file chosen

3. Select the 'Upload Data File' to upload the document.

- For the blank template, provide the following information: method code, other method name (if applicable), analyte code, run number (1,2,3,4,5), lot number, and results.
- For the pre-filled template, provide the following information: method code, other method name (if applicable), and results.

6. Do not alter column names, column order, or format of your template file. Any changes to the structure or format of the template will inactivate the template and the upload will fail.

7. Leave blank any inapplicable columns or fields. Add <LOD for results when necessary.

8. Upload the completed template for **ALL QC PROGRAMS** you are reporting. No partial results will be accepted at this time

9. **If you upload more than once, your previous file will be overwritten.**

10. Use the "Review QC Upload Data" page to review and edit entries (if necessary) after uploading the completed template.

11. If you have trouble uploading your completed template please email NSQAPDMT@cdc.gov for assistance.


Attach the Excel Template file *

Choose File
No file chosen


UPLOAD DATA FILE

4. After the file has been uploaded two notification emails will be sent:

a. Notification the QC file was uploaded.

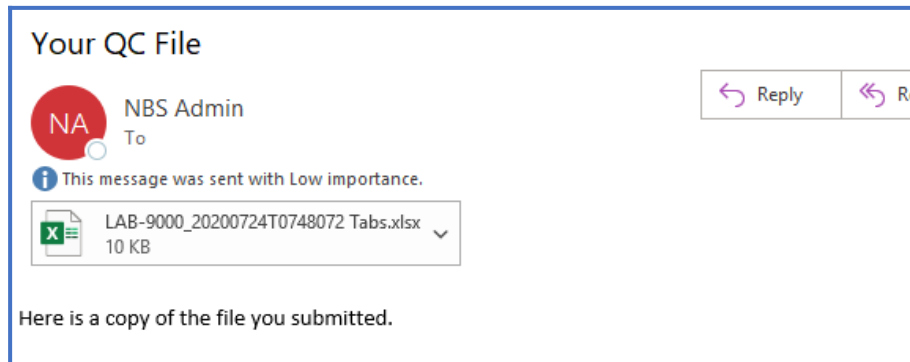


b. Notification data is ready for review in the Portal.

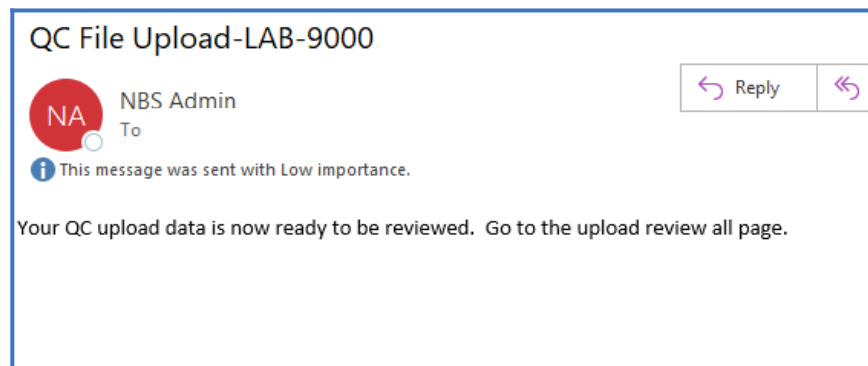


Note: It may take several minutes for this email to be sent. If your lab doesn't receive an email notifying the data is ready for review, an error with the upload may have occurred. Review Section 3.6 Troubleshooting and FAQs for more information.

3.3 Review Uploaded Data

1. Navigate back to the 'QC Data Upload Information' page and select the 'Review QC Upload Data' icon.

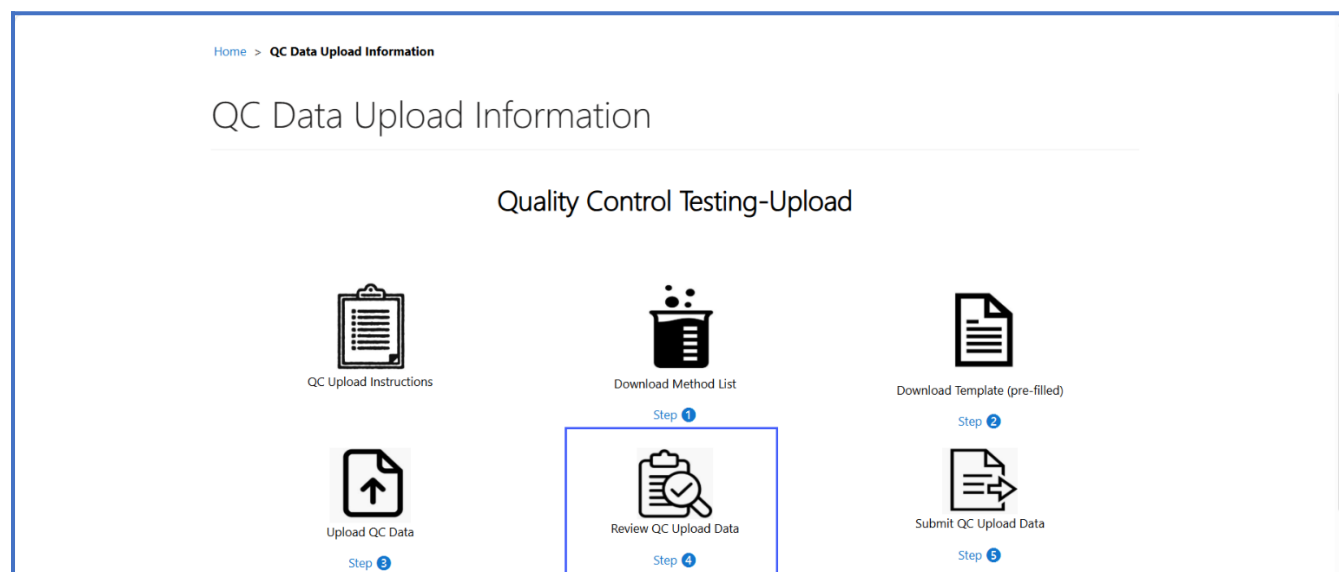

2. Review/edit the data that has been uploaded from the 'QC Review Upload Data' page.
 - a. Filter the data by program by selecting the check box next to the program(s) and selecting the 'Apply' button.

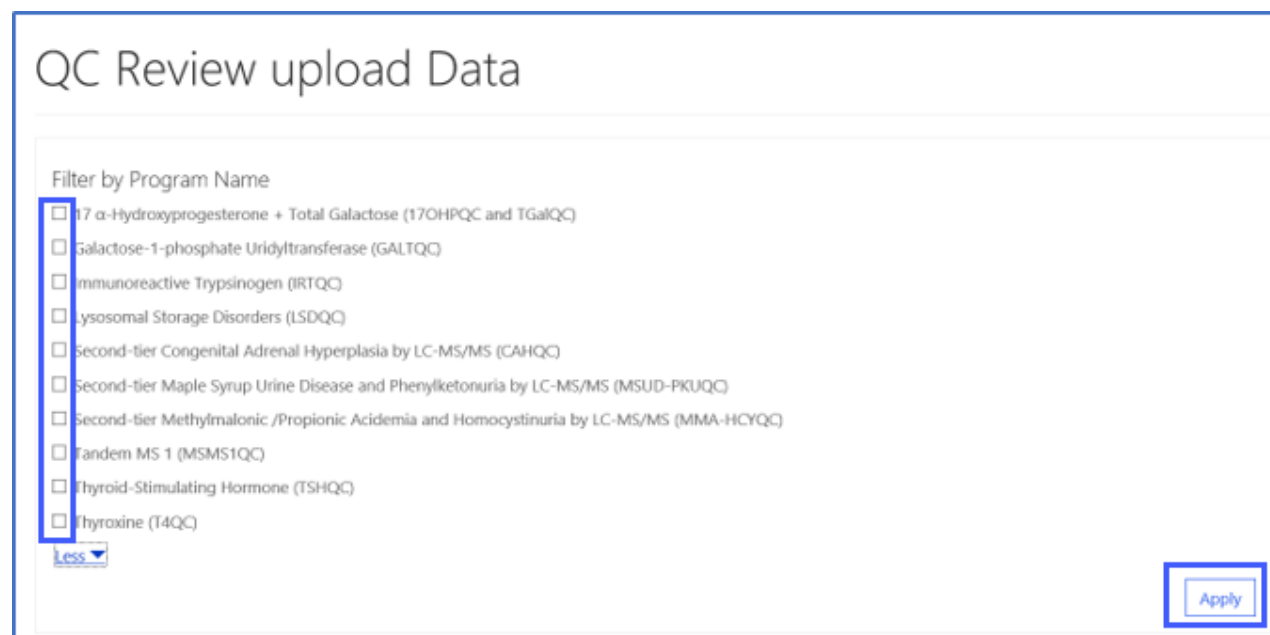

b. To sort data, click on the column header.

☐ Second-tier Congenital Adrenal Hyperplasia by LC-MS/MS (CAHQQ)
☐ Second-tier Maple Syrup Urine Disease and Phenylketonuria by LC-MS/MS (MSUD-PKUQQ)
[More](#)

[Apply](#)

[Download](#)

RUN	Method	Method Code	OTHER_METHOD	Analyte ↑	Analyte_Code	Analyte Abbreviation	Lot_A	Replicate_1A	Replicate_2A	Lot_B
1	LC-MS/MS non-kit	164		11-Deoxycortisol (11D2)	203	11D2	A1811	1.00	2.00	B1811
2	Other	999	DEMO	11-Deoxycortisol (11D2)	203	11D2	A1811		3.00	B1811
3	LC-MS/MS non-kit	164		11-Deoxycortisol	203	11D2	A1811			B1811

Note: The screenshot is sorted in ascending order by Analyte. The indicator is the arrow next to the header name.

c. To edit a record, select the number hyperlink in the 'Run' column.

☐ Second-tier Congenital Adrenal Hyperplasia by LC-MS/MS (CAHQQ)
☐ Second-tier Maple Syrup Urine Disease and Phenylketonuria by LC-MS/MS (MSUD-PKUQQ)
[More](#)

[Apply](#)

[Download](#)

RUN	Method	Method Code	OTHER_METHOD	Analyte ↑	Analyte_Code	Analyte Abbreviation	Lot_A	Replicate_1A	Replicate_2A	Lot_B
1	LC-MS/MS non-kit	164		11-Deoxycortisol (11D2)	203	11D2	A1811	1.00	2.00	B1811
2	Other	999	DEMO	11-Deoxycortisol (11D2)	203	11D2	A1811		3.00	B1811
3	LC-MS/MS non-kit	164		11-Deoxycortisol	203	11D2	A1811			B1811

- d. Edit the information on the pop-up screen for the run and select the **'Save Run Data'** button.

Analyte: 11-Deoxycortisol (11D2)

RUN 2

SAVE RUN DATA


Data cannot not be saved until all replicates have been completed

Values


A1811			
Replicate 1A	☒ <LOD	Replicate 2A *	☐ <LOD
—		3.00	
B1811			
Replicate 1B *	☐ <LOD	Replicate 2B *	☐ <LOD

3.4 Submit Uploaded Data


1. After data has been reviewed, navigate back to the 'QC Data Upload Information' page and select the '**Submit QC Data Upload**' icon.

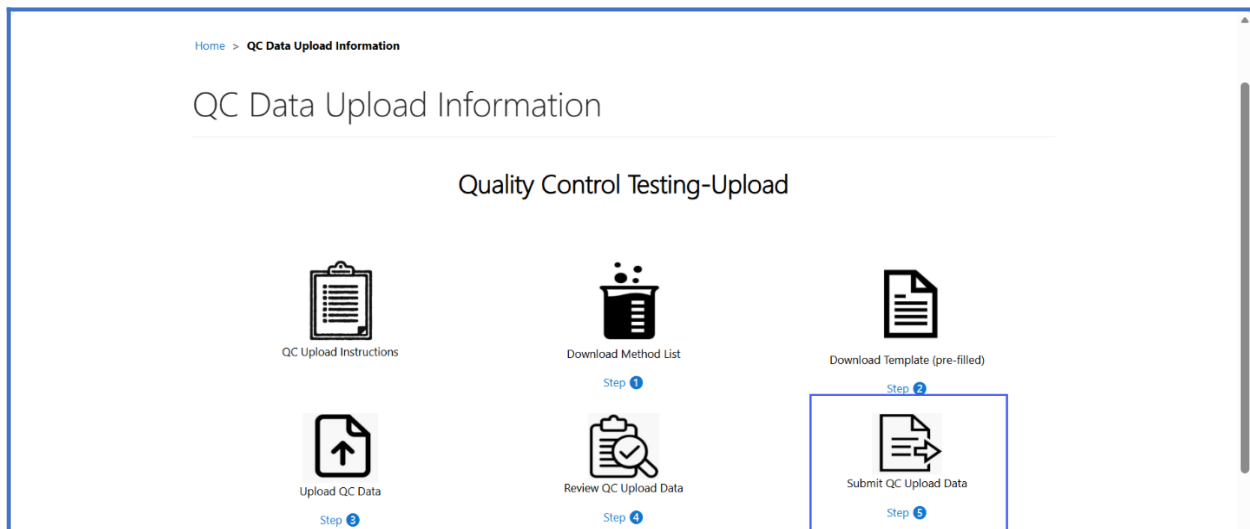

2. Select the '**Submit**' button to submit the uploaded data.

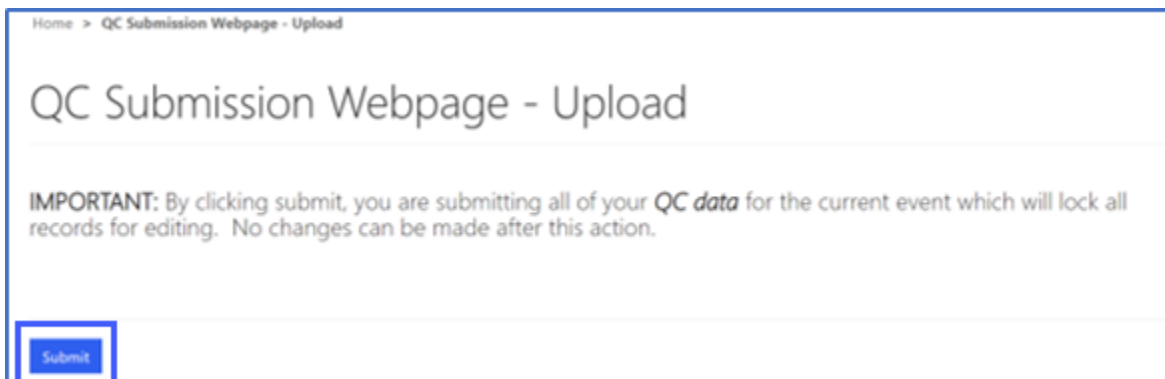

3. The submission confirmation message will appear upon successful upload.

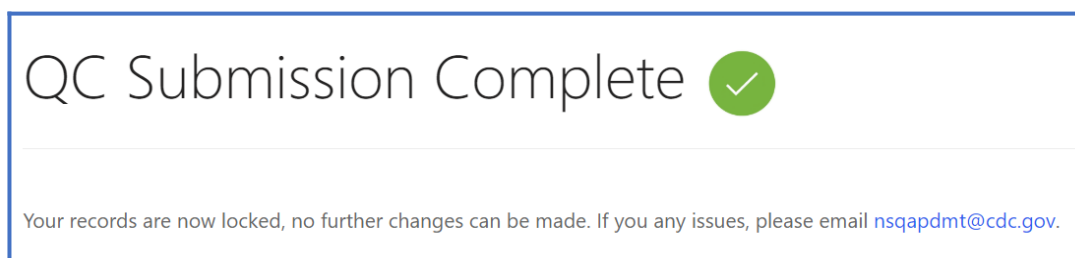

3.6 Troubleshooting and FAQs

If your uploaded data does not appear in the Portal after uploading

1. Confirm that you clicked the 'Review QC Upload Data' icon rather than the 'QC Manual Review' icon.
2. Confirm that you received a confirmation email stating that your data is ready. Note that it may take several minutes for this email to be sent.
3. If you received an email confirmation, refresh the web page.
4. If you clicked the 'Review QC Upload Data' icon, received an email confirmation, and refreshed the webpage and still do not see your data, confirm there are no issues with the file you uploaded. The following are common issues with the uploaded file:
 - a. The file should be in .xlsx or .xls format.
 - b. The data should appear on sheet 1 of the file.
 - c. Remove any table names in the spreadsheet. Table names are not encouraged unless you name the table 'Data'.
 - d. Check your headings. Make sure the headings are spelled correctly.
 - e. Check the values that have been entered for your data.
5. Confirm that your data was not already submitted.

If your method code does not appear

1. Confirm that the code exists in the method list that is provided on the Portal.
2. Check to see if any spaces or other characters are present before or after code. Remove as necessary.